

seca 656

Electronic platform
scales with innovative
memory function



- Especially spacious platform for all common roll-in systems.
- Stable construction with very high capacity.
- Transportable due to castors and integrated handles.
- Innovative memory function for determining precise weight of patient.

seca 656:

Precision weighing in all situations.

Emergency admission to hospital – everything has to happen very quickly. Especially when an operation is necessary and the patient’s weight must be determined in order to ensure the correct dose of anaesthetic. The weighing of patients on stretchers or gurneys in this or in less dramatic situations is no problem with the seca 656: the platform is large enough for all common gurneys or stretchers on undercarriages. And, of course, patients can also be weighed comfortably on a chair or in a wheelchair with this platform scale. The seca 656 is thus suitable for use in the dialysis sector and in old people’s and nursing homes.



Stable construction with very high capacity.

Thanks to its extremely stable construction and very high capacity of 360 kg, even heavyweight patients can be weighed easily with the platform scale 656. The platform is generously sized to cope with wheelchairs and all common roll-in systems and stretchers on undercarriages.

Castors and integrated handles for transport.


The seca 656 is a very versatile scale as it is easy to move around on its smooth-running castors with the aid of the integrated handles.

Cable remote display with many additional functions.

Included in delivery is the cable remote display with easy to read LCD figures. It can be mounted on a table or hung on a wall at eye level. The display can be operated intuitively and is equipped with many intelligent functions: Thanks to the HOLD function, the patient can first be attended to before the result is noted. The innovative memory function of the seca 656 is particularly practical – for example in the casualty department before an operation: the total weight of the patient including stretcher or

wheelchair is first determined. After again weighing the stretcher or wheelchair the net weight of the patient appears automatically. With the pre-TARE function, the additional weight of auxiliary aids can also be deducted from the total weight while actually weighing the patient. And, finally, the individual nutritional condition of the patient can be ascertained with the aid of the BMI function.

With its transport castors and steering castor, the seca 656 can be maneuvered easily.

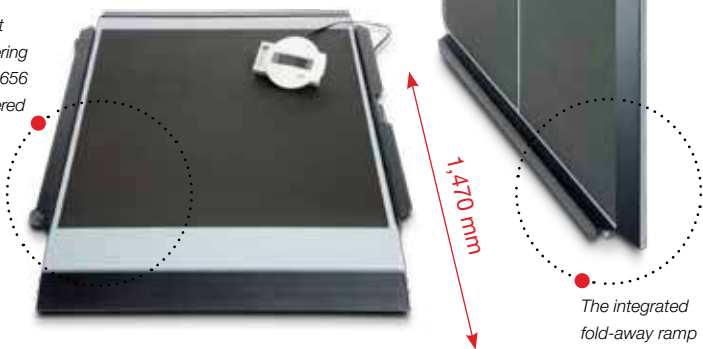

The integrated fold-away ramp locks into place for troublefree transport.

seca 656

Technical Data

- Capacity: 360 kg / 800 lbs
- Graduation: 100 g / 0.1 lbs
- TARE-range: 360 kg / 800 lbs
- Dimensions (WxHxD): 893 x 87 x 1,655 mm / 35.2 x 3.4 x 65.2"
- Dimensions / platform (WxHxD): 760 x 55 x 1,470 mm / 29.9 x 2.2 x 57.9"
- Weight: 46.2 kg / 101.9 lbs
- Power supply: Power adapter
- Functions: TARE, pre-TARE, auto-HOLD, BMI, kg/lbs switch-over, damping
- Optional: stand for cable remote display seca 472, RS232 interface seca 460