

seca 869

Flat scales with cabled remote display



- Separate display and operating element with infinitely adjustable tilt angle for convenient use.
- Large, easy-to-read 3/4-inch-high LCD figures.
- BMI function to determine nutritional status.

seca 869:

What was difficult is now easy.

The seca 869 is exactly right for those looking for a scale with practical cabled remote display and a handle for flexible use. This scale has many functions that make life easier. It's simple to use and easy to clean.



Read the display while at your desk.

The seca 869 makes your work easy. Simply place the display on your desk and the scale in front of it. Connect the two with the 7¼ ft. cable and you're ready to work. It's also easy to read the results since the tilt angle of the display is infinitely adjustable.

Stability even on problematic surfaces.


An especially large leveling base gives the seca 869 great stability. The display and operating element is similarly equipped with a stable rubber base and with an easy-to-use self-explanatory keyboard. Storage is also effortless, thanks to its compact dimensions and light weight.

Measuring with all essential functions.

The seca 869 has several convenient functions required in routine medical practice: BMI to determine nutritional status, HOLD to keep the measurement results in the display, TARE to tare the weight, e.g., to weigh a small child held in a parent's arms.



Leveling base on the bottom of the scale base for safe and secure placement.

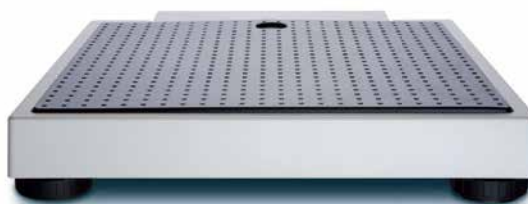

Compact design and low weight for great carrying comfort and universal usage.



Infinitely adjustable tilt angle in display and operating element for installation on a desk or wall.



seca 869

Technical Data

- Capacity: 250 kg / 550 lbs
- Graduation: 100 g < 150 kg > 200 g / 0.2 lbs < 330 lbs > 0.5 lbs
- Dimensions (WxHxD): 321 x 60 x 356 mm / 12.6 x 2.4 x 14"
- Dimensions/platform (WxHxD): 280 x 60 x 278 mm / 11 x 2.4 x 10.9"
- Cable length: 2 m / 78.7"
- Weight: 4.5 kg / 9.9 lbs
- Power supply: Batteries, power adapter optional
- Functions: TARE, HOLD, BMI, lbs/kg switch-over, automatic switch-off
- Optional: stadiometer seca 217 combined with adapter element seca 437, switch-mode power adapter seca 400, carrying case seca 414 (for combination with stadiometer and adapter element), carrying case seca 421