


SPECTRO₂ | LOGIX™

Clinical Analysis Software

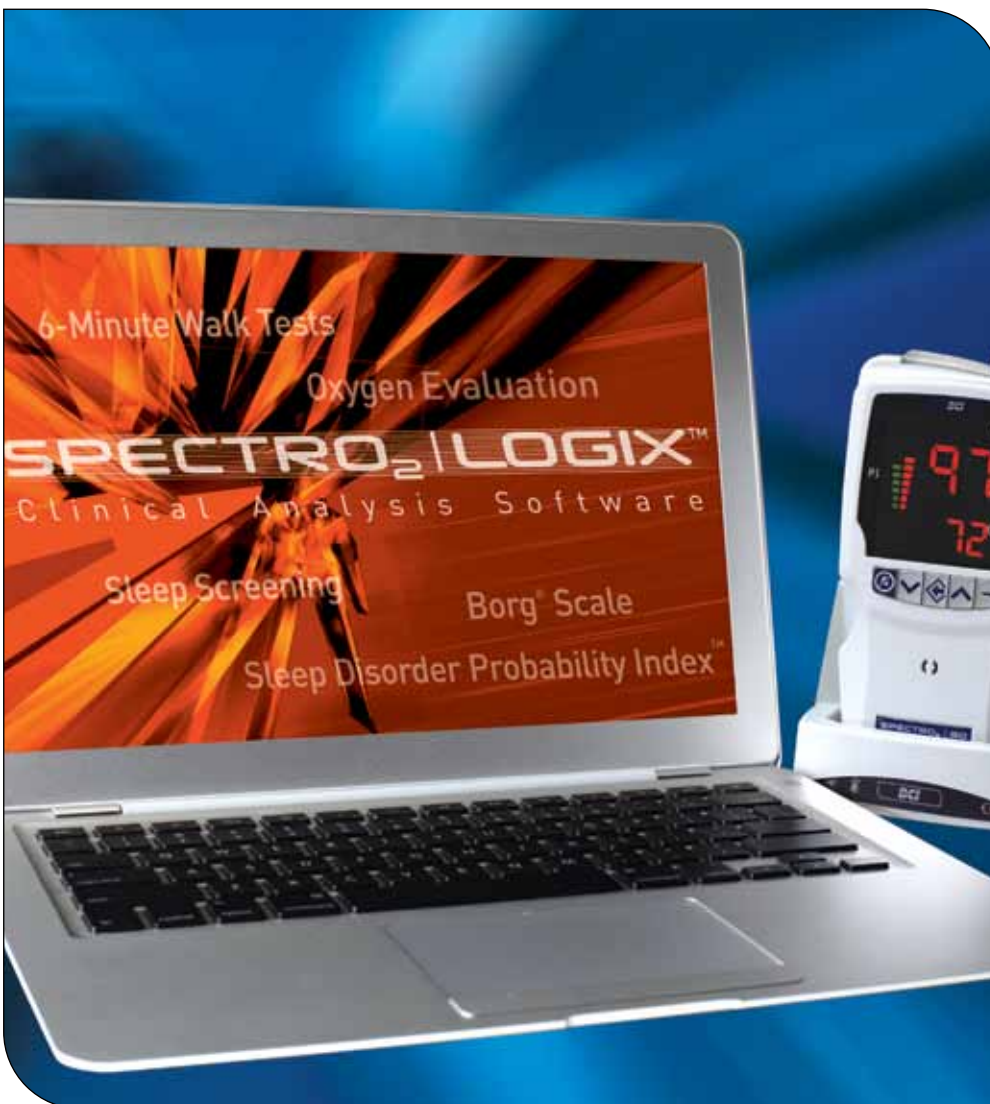

SPECTRO₂ | LOGIX™ Clinical Analysis Software streamlines pulse oximetry data collection and helps eliminate documentation errors during:

- Sleep screening
- Oxygen evaluation
- 6-minute walk tests
- Ambulatory analysis

Designed with clinicians for clinicians, SPECTRO₂ | LOGIX™ Clinical Analysis Software provides more diagnostic analysis capabilities than other pulse oximetry software, including:

- The proprietary Sleep Disorder Probability Index, or SDPI™ assessment tool – to indicate the relative potential a patient may suffer from sleep disordered breathing
- Borg® scale – helps to measure a patient's perceived rate of exertion without manual calculations

PATIENT MONITORING

**SPECTRO₂ | LOGIX™ software
designed with you in mind**

SPECTRO₂ | LOGIX™

Clinical Analysis Software

Ambulatory and Six Minute Walk Tests (6MWT)


The **Ambulatory** and **Six Minute Walk Tests (6MWT)** provide the clinician with the ability to record oxygen saturation and heart rate levels to assess the patient's oxygen needs.

The **6MWT** is a very simple diagnostic test to help assess a patient's exercise tolerance over a period of six minutes. The patient walks at his own intensity as the data is captured by the SPECTRO₂™ pulse oximeter. The data is then transferred to SPECTRO₂ | LOGIX™ to quickly help the clinician assess the relative degree of exertion as perceived by the patient by utilizing a modified **Borg® Scale**.

Through innovation and engineering superiority in the field of Pulse Oximetry and Sleep Screening, the BCI® brand of Smiths Medical offers SPECTRO₂™ Pulse Oximeters, a comprehensive Pulse Oximetry portfolio designed to assist clinicians with patient monitoring ranging from day-to-day wellness checks to more in-depth patient monitoring for efficacy of respiratory treatments.

The SPECTRO₂ | LOGIX™ Clinical Analysis Software application is a unique pulse oximetry software tool designed for exclusive use with the SPECTRO₂™ Pulse Oximeters. Serving as a centralized data management and analysis system, it captures demographic, visit, medication history and pulse oximetry testing information for each patient, seamlessly delivering analysis and reporting capabilities of the pulse oximetry data. Customizable reports allow the clinician to select, zoom in on, and even remove specific data segments of the analysis. After identifying and capturing the pertinent test information, a report can be printed and used for further analysis by other clinicians or submitted for reimbursement claims.

Features such as a modified **Borg® scale** and the **Sleep Disorder Probability Index, or SDPI™ assessment tool**, which is proprietary to the BCI® SPECTRO₂ | LOGIX™ software, help to accelerate patient assessments by providing immediate analysis to assist in determining the patient care plan in a timely manner. The automatic transfer of pulse oximetry data to SPECTRO₂ | LOGIX™ software can also help reduce paperwork and documentation errors, and the analysis data can be stored for future use.

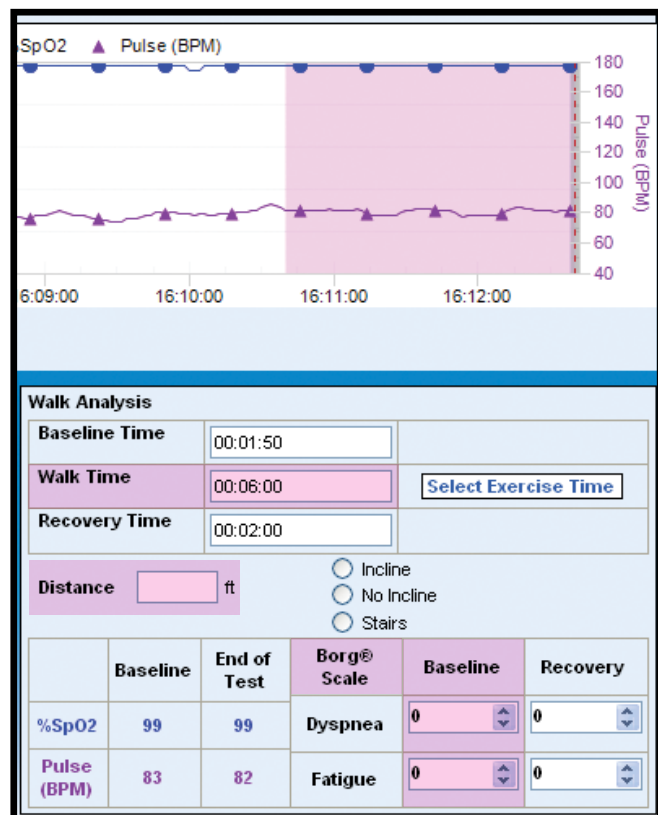

Enter other pertinent information such as baseline readings, time and distance walked.

Sleep Screening and Oxygen Evaluation

The **Oxygen Evaluation Test** is usually performed in the patient's home during sleep. A clinician will show the patient how to use the pulse oximeter. Once captured, the data can be downloaded to the application on a PC in the physician's or testing facility's offices for evaluation.


Duration of the longest occurrence with %SpO2 <= 88:		00:00:00					
Start time for this occurrence:		--					
Total Number of Events:		0					
Event	Start Date / Time	Stop Date / Time	Duration	Min %SpO2	Avg. BPM	Max. BPM	Min. BPM

This analysis shows Mean and Longest Durations of desaturation events. It also lists the 50 deepest desaturation events, their minimum %SpO₂, minimum and maximum pulse rate and duration.

Sleep Screening is used to screen for sleep disordered breathing, primarily Obstructive Sleep Apnea (OSA).

Based on criteria identified by the American Academy of Sleep Medicine, SPECTRO₂ | LOGIX™ Clinical Analysis Software utilizes its proprietary **Sleep Disorder Probability Index, or SDPI™ assessment tool**, to assist in determining a patient's probability of having a Sleep Disorder.

The **SDPI™ assessment tool** rates the relative probability that the patient may suffer from sleep-disordered breathing and require further testing in the sleep lab.

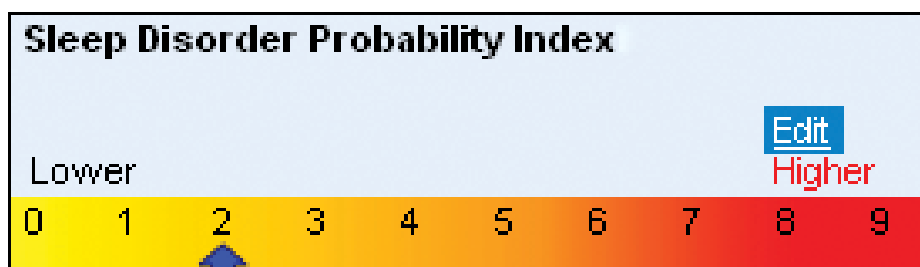

Rates the relative probability that the patient may suffer from sleep-disordered breathing and require further testing in the sleep lab.

The SDPI™ assessment tool is one of the most unique features of the SPECTRO₂ | LOGIX™ Clinical Analysis Software application.

The BCI® SPECTRO₂™ Pulse Oximeters

SPECTRO₂ | LOGIX™ Clinical Analysis Software is for exclusive use with the Smiths Medical BCI® brand of SPECTRO₂™ Pulse Oximeters.

- **SPECTRO₂™ | 10 pulse oximeter**
Reliable results at an affordable cost.
- **SPECTRO₂™ | 20 pulse oximeter**
Patented BCI® serial autocorrelation technology delivering reliable results in challenging clinical settings such as low perfusion and motion.
- **SPECTRO₂™ | 30 pulse oximeter**
Patented BCI® brand serial autocorrelation technology with low perfusion performance and motion tolerance. Audible and visual alarms, as well as a remote nurse call alarm allow healthcare providers to monitor remotely and help promote patient safety.

When monitored real-time by a clinician, the **SPECTRO₂™ | 10** and **SPECTRO₂™ | 20** pulse oximeters are ideal for the **6MWT** and other **ambulatory analyses**.



SPECTRO₂™ | 30 Pulse Oximeter is used to monitor the efficacy of respiratory treatments, and is ideal for use in obtaining **Oxygen Evaluation** and **Sleep Screening analyses**.

For additional information on the SPECTRO₂™ Pulse Oximetry portfolio, please visit www.smiths-medical.com

SPECTRO₂ | LOGIX™ Clinical Analysis Software System Requirements:

The following is required to run the software:

- Pentium® Processor
- Windows® XP Pro (Service Pack 3) 32-bit (x86)
- .NET Framework 2.0 (automatically installs when software installs)
- Crystal Reports for .NET Framework 2.0 (automatically installs when software installs)
- Microsoft® SQL Server® 2005 Express Edition SP2 (automatically installs when software installs)
- 1 GB RAM minimum
- 1 GB hard drive space minimum
- Screen resolution of 1024 x 768 required (1280 x 1024 resolution recommended)

Catalog Number – 21-0469-0105-24 Includes:

- 21-0470-0105-24 SPECTRO₂ | LOGIX™ Clinical Analysis Software (CD ROM)
- WW1089 USB Cable

ECOMED

Serving The HealthCare Industry

t: 1800 675 432

f: 1800 675 433

e: enquiries@ecomед.com.au

w: www.ecomed.com.au

Smiths Medical ASD, Inc.
Norwell, MA 02061, USA
Phone: 1-214-618-0218
Toll-Free USA: 1-800-258-5361
www.smiths-medical.com

Find your local contact information at:
www.smiths-medical.com/customer-support

BCI, SPECTRO₂, SPECTRO₂ | LOGIX, SDPI and the BCI and Smiths Medical design marks are trademarks of Smiths Medical. The symbol ® indicates the trademark is registered in the U.S. Patent and Trademark Office and certain other countries. All other names and marks mentioned are the trade names, trademarks or service marks of their respective owners. ©2011 Smiths Medical. All rights reserved. PM195636GL 09/11

smiths medical