

ULTRAVIEW SL

SvO₂ / ScvO₂ MODULE (91424)



The Ultraview SL™ SvO₂ / ScvO₂ module measures mixed venous oxygen saturation (SvO₂) or central venous oxygen saturation (ScvO₂), providing a status indicator of oxygen delivery and consumption in critically ill patients (see Figure 1).

Clinical Benefits

- Measurement of ScvO₂ with a central venous catheter, rather than SvO₂ with a pulmonary artery catheter, may allow a broader range of patients to be monitored.
- SvO₂ / ScvO₂ values and trends provide clinicians with information to assess cardiovascular and respiratory compromise.
- May be used as a component of sepsis screening and early goal-directed therapy protocols.
- SvO₂ / ScvO₂ saturation and light intensity values are continuously added to the patient record, making it easy to screen critical patients using the ICS G2 Custom Trends feature (see Figure 2).

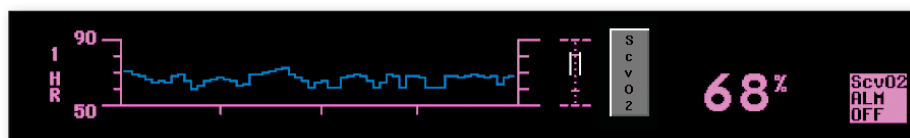


Figure 1



Figure 2


Ultraview SvO₂ / ScvO₂ Module (91424)

Features

- Provides continuous, real-time monitoring using proven TriOx™ 3 wavelength oximetry from ICU Medical.
- Modular design allows clinicians to perform SvO₂/ScvO₂ monitoring at any bed as needed.
- Provides oximetric calculations for:
 - CaO₂
 - CvO₂
 - O₂AV
 - avDO₂
 - VO₂

Supplies and Accessories

The Ultraview SL SvO₂ / ScvO₂ module is compatible with ICU Medical's TriOx Optical Module, fiber-optic SvO₂ / ScvO₂ catheters and ScvO₂ probes*.

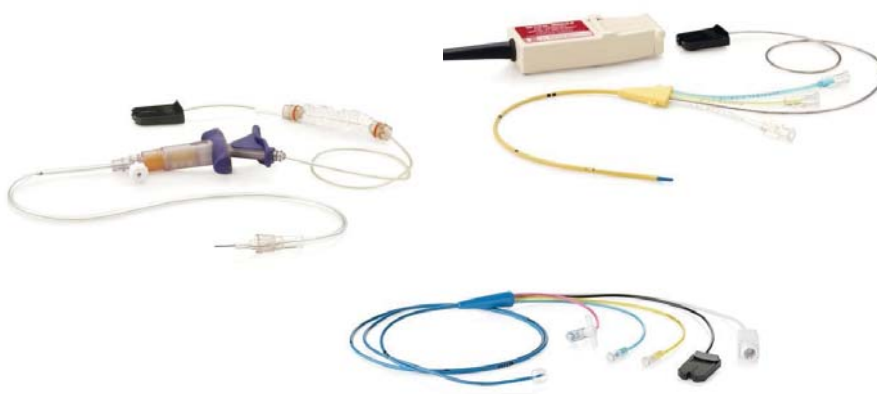

ICU Medical | 1.800.824.7890 | www.icumed.com

* May eliminate the need to exchange central venous catheters in order to monitor ScvO₂ continuously. The ScvO₂ probe can be inserted into an existing 16 cm or 20 cm central line with a minimum distal lumen inner diameter of 0.032" (16 G) and a maximum length from the outside of the female luer on the distal lumen to the catheter tip of 35 cm. Contact ICU Medical for more information.

Distributed by

ECOMED

1800 67 5432

1800 67 5433

sales@ecomед.com.au

www.ecomed.com.au

www.spacelabshealthcare.com



Spacelabs Healthcare Company Headquarters
5150 220th Ave. SE, Issaquah, WA 98029
PO Box 7018, Issaquah, WA 98027-7018

Telephone: +1 425-657-7200
In North America: +1 800-522-7025
Fax: +1 425-657-7212

© 2012 Spacelabs Healthcare
Specifications subject to change without notice. Ultraview SL is a trademark of Spacelabs Healthcare. TriOx is a trademark of ICU Medical. Other brands and product names are trademarks of their respective owners.

Product may not be available in all locations. Check with your local Spacelabs representative.

M1016-my12