


Data you can rely on.
People you can trust.



Vitalograph Pneumotrac with Spirotrac Software

Combines the power of your Windows® PC with Vitalograph measuring technology for early detection of respiratory disease and integration with your institution network.

- Accurate, robust and linear Fleisch pneumotachograph
- Intuitive user interface with workflow orientated presentation
- Automatic FEV₁ trend chart of actual and normal range data as subject is selected
- Open session, f/v and v/t curves, all test data and test quality in one view
- Automatic storage of all test data and events, including calibration log
- Built-in test quality checks
- Choice of exciting incentives
- Fully scalable database capacity using powerful Microsoft® SQL server technology
- Network ready

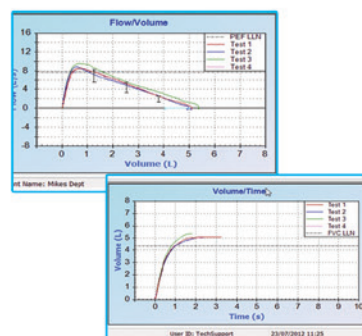

Vitalograph Spirotrac

Cardio-respiratory diagnostic software

- Network multiple Spirotrac® cardio-respiratory diagnostic software devices in their own network or to a central database on your existing networks
- Full spirometry suite
- 12-lead ECG Bluetooth®
- Pulse Oximetry USB
- Bronchial Challenge Testing
- 6 Minute Walk Test and COPD Assessment
- Blood Pressure USB
- Weighing Scale Bluetooth



Technical Specification:

Product: Pneumotrac™
Model No: 6800
Size: 170mm x 163mm x 45mm
Net Size: 103mm x 185mm
Weight: Net 450g Gross 650g
Operating Temperature Range: ATS/ERS limits: 17-37°C Design Limits: 10-40°C
Power Supply: From USB 2 port
Flow Detection Principle: Fleisch type pneumotachograph
Volume Detection: Flow integration sampling @ 100Hz
Accuracy when operated in operating temperature range conditions: Volumes: Better than ±3% (Max displayed volume 10L)
Flows: Better than ±10%(Max ±16 L/s / Min ± 0.02 L/s)
Max Test Duration: 45 seconds
Environmental Data: Selectable temperature (built in sensor), barometric pressure, altitude & humidity
Parameters Displayed: User selectable from: VC, IVC, IC, VT (TV), TLC, RV, IRV, ERV, FRC, FVC, FIVC, FEV₁, FEV₃, FEV₆, FVC, FEV₁/VC, FEV₁/FVC, FEV₃/VC, FEV₃/FVC, FEV₁/FEV₆, FEF₇₅, FEF₅₀, FEF₂₅, FEF₂₅₋₇₅, FEF₂₅₋₇₅/FVC, FIV₁, PIF, FIV₁/FIVC, FIF₂₅, FIF₅₀, FIF₇₅, FEF₅₀/FIF₅₀, FET, MVVind, FEV₁ Ratio, FEV₀₋₅, PEF L/min, PEF L/s, FEF₀₋₂₋₁₋₂, FEF₇₅₋₈₅%, FEF₂₅%, FEF₅₀%, FEF₇₅%, FMFT, FIF₂₅%, FIF₅₀%, FIF₇₅%, PIF, FEV₀₋₇₅, and lung age
Test Types: Single breath tests, flow/volume loops, multi-breath testing, tidal breathing and combined VC/FVC type test methods supported
Predicted Values Selectable: ERS/ECCS, NHANES, Polgar, Pereira, Berglund, Forsche, Gutierrez, Hedenström, Taiwan, Knudson, Crapo, Hsu, KNLW, Viljanen, SEPAR, Gulsvik, SBPT, Tamura, Ip, Quanjer, Wang, Dockery and more
Performance Standards: EN ISO 23747:2009 & ISO 26782:2009 & ATS/ERS2005
Safety Standards: EN 60601 (IEC 60601)
Medical Safety Standard: Medical Devices Directive 93/42/EEC-1993
QA/GMP Standards: ISO 13485:2003, FDA 21 CFR 820, CMDR SOR 98/282 and JPAL

Ordering Info:

77700 6800 Pneumotrac with Spirotrac®
77550 6800 Pneumotrac with Spirotrac & ECG
70201 7000 Spirotrac Network Licence

Note: Spirotrac Software includes software for Spirometry, 12 Lead ECG, Pulse Oximetry, 6MWT, NIBP, Weight and Bronchial Challenge Testing

Compatible Devices and Hardware:

40400 4000 Respiratory Monitor asma-1™ USB
41300 12 Lead Bluetooth ECG
70400 Pulse Oximetry Finger Probe USB
70401 Biometric Scanner USB
70402 Blood Pressure Monitor USB
70404 Medical Weighing Scales Bluetooth

Consumables and Accessories:

41301 ECG Electrodes Disposable (30)
28350 2820 BVF™ Bacterial/Viral Filters (50)
20242 2024 SafeTway® Mouthpieces (200)
20980 2024 SafeTway® Mini Mouthpieces (50)
20303 2030 Noseclip Disposable (200)
36020 2040 Precision Syringe 3-L
42084 Flow Conditioning Meshes (10)

Vitalograph is a world leading provider of outstanding quality cardio-respiratory diagnostic devices, clinical trial services and medical equipment servicing. With a pioneering heritage of excellence spanning half a century Vitalograph continues to make valuable contributions to effective medical care and enhanced quality of life.

Distributed by



1800 67 5432



1800 67 5433



sales@ecomед.com.au



www.ecomed.com.au



Your cardio-respiratory partner

Vitalograph Ltd. UK
Maids Moreton, Buckingham,
MK18 1SW, England
Phone: (01280) 827110
Fax: (01280) 823302
e-mail: sales@vitalograph.co.uk
www.vitalograph.co.uk

Vitalograph Ltd. Export
Maids Moreton, Buckingham,
MK18 1SW, England
Phone: +44 1280 827120
Fax: +44 1280 823302
e-mail: sales@vitalograph.co.uk
www.vitalograph.eu

Vitalograph GmbH
Rellinger Straße 64a,
D-20257 Hamburg,
Germany
Phone: (040) 547391-0
Fax: (040) 547391-40
e-mail: info@vitalograph.de
www.vitalograph.de

Vitalograph Inc.
13310 West 99th Street,
Lenexa, Kansas, 66215,
USA
Toll Free: 800 255 6626
Phone: (913) 730 3200
Fax: (913) 730 3232
e-mail: vitcs@vitalograph.com
www.vitalograph.com

Vitalograph (Ireland) Ltd.
Gort Road Business Park,
Ennis, Co. Clare,
Ireland
Phone: +353 65 6864100
Fax: +353 65 6829289
e-mail: sales@vitalograph.ie
www.vitalograph.ie

Vitalograph Ltd. Hong Kong
Unit 230, Corporation Park,
11 On Lai Street,
Shatin,
Hong Kong
Phone: +852 2117 2678
Fax: +852 2117 2679
e-mail: sales@vitalograph.hk
www.vitalograph.hk